

AB-PINACA N-Pentanoic Acid ELISA Standard

Cat. #9028 LotQ0815


LIMITATIONS: THIS PRODUCT IS FOR RESEARCH USE ONLY AND IS NOT APPROVED FOR THERAPEUTIC OR DIAGNOSTIC USE.

Applications:

The AB-PINACA N-pentanoic acid ELISA Standard is intended to be used for preparing a control standard for Cat. #4800 AB-PINACA Synthetic Cannabinoid ELISA Kit. The product must be diluted in an appropriate buffer to the desired concentration for use.

Product Description:

AB-PINACA N-pentanoic acid: 5-(3-((1-amino-3-methyl-1-oxobutan-2-yl)carbamoyl)-1H-indazol-1-yl)pentanoic acid
CAS Number: N/A
Formula Weight: 360.4
Structure:

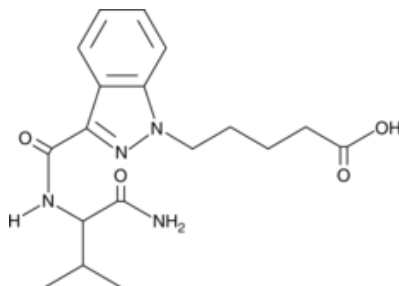

Supplied As:

0.5mL of 1 microgram/mL AB-PINACA N-pentanoic acid in phosphate buffered saline (20mM sodium phosphate, 150mM NaCl, pH7.4)

Storage and Stability:

Stable for 6 months from date of shipment when stored at -20°C (avoid frost-free freezers) or >1year at -70°C.

Usage:

- 1) Defrost #9028 product tube at room temperature.
- 2) Vortex tube vigorously for 20 sec.
- 3) Pipette an aliquot of #9028 into the desired buffer or sample matrix (e.g.: phosphate buffered saline, synthetic urine).
- 4) Vortex vigorously for 20 sec.
- 5) Use diluted sample immediately or aliquot and store frozen at -20 or -70°C.

For example, to make a 10ng/mL dilution:
9.9mL PBS or sample matrix
100µL AB-PINACA N-pentanoic acid,
1µg/mL

Tulip BioLabs Related Products:

- Cat. #4800 **AB-PINACA Synthetic Cannabinoid ELISA Kit**
- Cat. #9023 **JWH-018 N-Pentanoic Acid ELISA Standard**
- Cat. #9024 **JWH-250 N-Pentanoic Acid ELISA Standard**
- Cat. #9025 **UR-144 N-Pentanoic Acid ELISA Standard**
- Cat. #9026 **PB-22 N-Pentanoic Acid ELISA Standard**
- Cat. #9027 **AKB48 N-Pentanoic Acid ELISA Standard**